

Your Guide for a Successful First Fit.



A strong first fitting builds trust with your clients from the very first sound. Phonak supports this moment with **clinically driven solutions and an easy-to-navigate fitting workflow**. Follow these **5 key steps** for a fitting with Target to make that first impression count.

Step 1

Connect the devices, allocate sides and select experience level.

Optimizing the outcome:

- Identify the device and determine which is Right/Left through the music icon or by pushing on the rocker switch.
- Initial EON fittings may require a manual Fitting authorization. If needed it only takes 8 seconds, and you will be guided through the process.

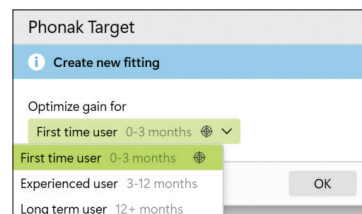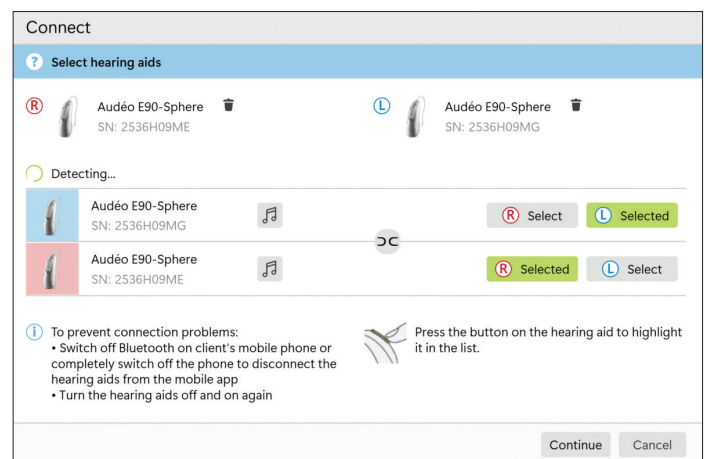


Step 2

Ensure the earpiece (dome or custom earpiece) in Target matches what is on the device.

Optimizing the outcome:

- The acoustics selected here will determine the intensity and duration of the feedback test signal as well as other feature settings.

Note: When using NOAH, the audiogram will automatically be pulled into Target. In Target standalone, it can be entered into the Audiogram screen.

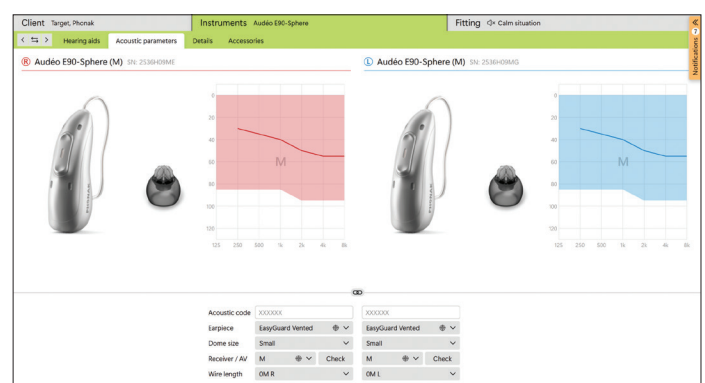

Step 3


Run the Feedback & real ear test to provide accurate information regarding how the earpiece is functioning in the client's ear. The amount of stable gain for this client can then be easily visualized.

Optimizing the amount of stable gain:

- The feedback & real ear test shows how much stable gain is available without activating the feedback algorithm.
- If sufficient gain cannot be provided, you have two options:
 1. Change to a more occluding earpiece.
 2. Apply more Additional headroom.

Note: EON devices will automatically apply up to 6 dB Additional headroom when required.

- To maintain sound quality, only apply enough Additional headroom so that the gain targets can be reached.

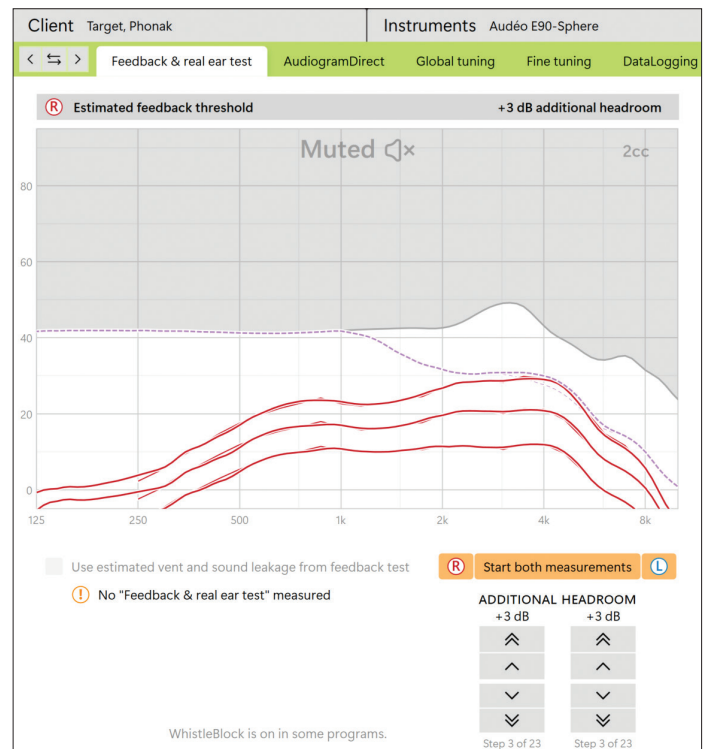


Step 4

Check the amplification and MPO with your client.

Optimizing the outcome:

- Verifying the fitting with Real Ear Measures allows you to provide a completely personalized fitting.
- If Auto Acclimatization is active, we recommend setting the master gain to 100% prior to making any adjustments to the Gain & MPO settings.
- The master gain can be set in Global tuning to meet your client's current needs. As your client acclimatizes to the amplification, this setting can be adjusted. The master gain setting influences the gain in all programs.

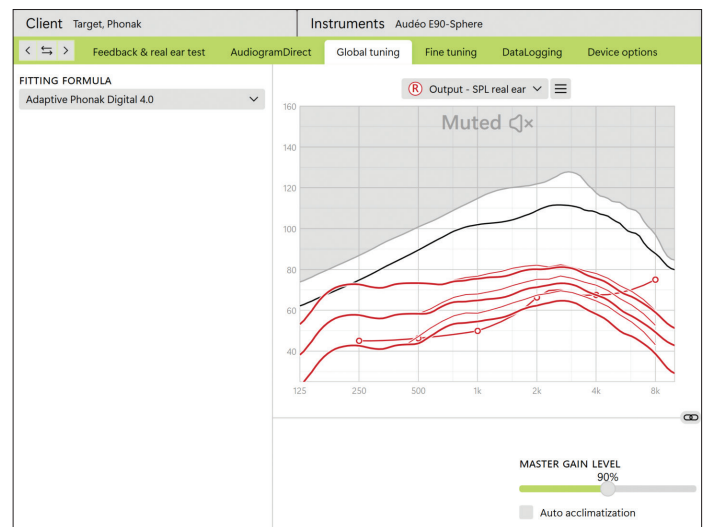

Step 5

Save the fitting to the hearing aids and Target session.