Roger On Sphere™

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Charge! | Laden! | Chargez! | Laden!

Importante Information

For Roger On to function properly your hearing device must be fitted with a Roger receiver. If the receiver is not attached/ installed, please refer to the enclosed quick guide "Roger receiver setup". If you are still unsure how to setup the receiver, we recommend contacting your hearing care professional for assistance.

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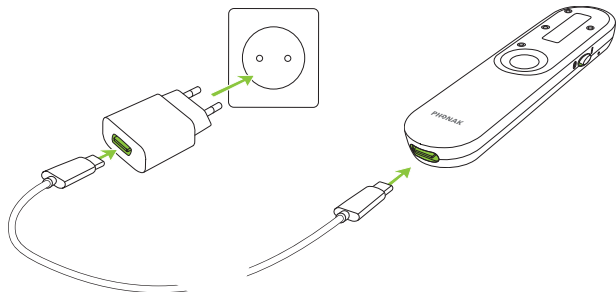
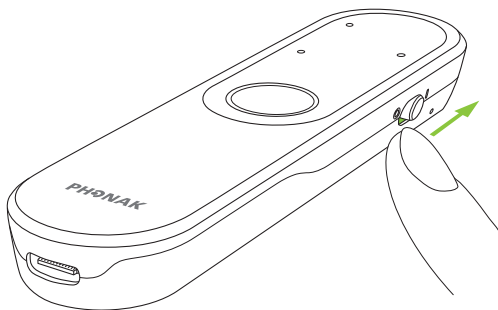
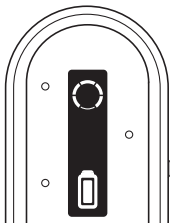
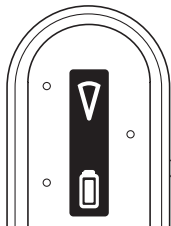
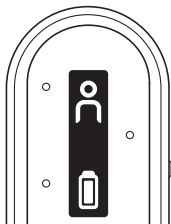
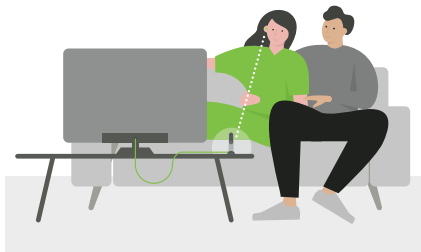
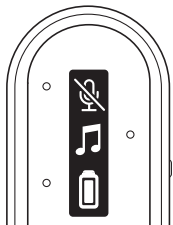
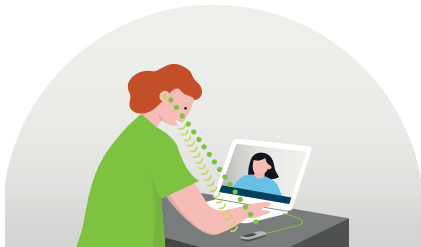
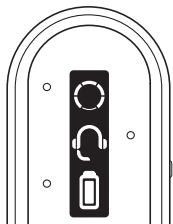
myRogerMic App

The Apple App Store logo is a trademark of Apple Inc.
The Google Play store logo is a trademark of Google LLC.

For further information, please refer to the User Guide. | Weitere Informationen zur Anwendung finden Sie in der Gebrauchsanweisung. | Pour plus d'informations, veuillez vous référer aux instructions d'utilisation. | In de gebruiksaanwijzing vindt u meer gedetailleerde informatie.

Compliance information

Declaration of Conformity

Hereby Phonak Communications AG declares that this product is in compliance with the Radio Equipment Directive 2014/53/EU.

The full text of the EU

Declaration of Conformity can be obtained from the manufacturer or the

local Phonak representative or at www.phonak.com

Operation frequency: 2400 – 2483.5 MHz

Power level: < 100 mW

Notice 1:

This device complies with Part 15 of the FCC Rules and with Industry Canada's license exempt RSS standard(s). Operation is subject to the following two conditions:

1. this device may not cause interference, and
2. this device must accept any interference, including interference that may cause undesired operation of the device.

Notice 2:

The device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. The device generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the device does cause harmful interference to radio

or television reception, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Notice 3:

Changes or modifications made to this device not expressly approved by Phonak Communications AG may void the FCC authorization to operate this device.

Notice 4:

FCC / Industry Canada RF Radiation Exposure Statement. This device complies with FCC RF radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

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WARNING

- Using this device may reduce background noises. Please be aware that warning signals or noises, e.g. cars, are partially or entirely suppressed.
- Do not use your device in areas where electronic equipment is prohibited.
- Do not use the device in explosive areas (mines or industrial areas with danger of explosions, oxygenrich environments or areas where flammable anesthetics are handled).
- The system operates between 2.4 GHz and 2.48 GHz frequency range. When flying please check if flight operator requires devices to be off.
- Protect the device from excessive moisture (bathing or swimming areas) and heat sources (radiator). Protect the device from excessive shock and vibration.
- If the device has been dropped or damaged, if it overheats, has a damaged cord or plug, or has been dropped into liquid, stop using your device and contact authorized service centre.
- Always make sure that the device and docking station are dry and clean.
- Clean the device using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean them. never use a microwave or other heating devices to dry the device.
- Store the device in a dry place if you are not using it for a long period of time.
- Changes or modifications to any of the devices that were not explicitly approved by Sonova Communications AG are not permitted. Such changes may damage your ear or the devices.
- Only authorized personnel are allowed to open the product and its accessories due to possible electrical hazard.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the earbuds or charging case, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
- This device may generate magnetic fields. If any influence of the device on the implantable device (e.g. pacemakers, defibrillators, etc.) is experienced, stop using the device and contact your doctor and / or the manufacturer of the implantable device for advice.
- Don't expose your device to temperatures outside the operation conditions indicated in the User Guide (see chapter Specifications). They contain lithium-ion batteries, which might swell, leak or explode upon exposure to temperatures outside those conditions.
- This device uses a lithium battery. These batteries are hazardous and can cause severe or fatal injuries in 2 hours or less if they are swallowed or placed inside any part of the body, whether the battery is used or new! Keep out of the reach of children, individuals with cognitive impairment, or pets. If you suspect that a battery was swallowed or placed inside any part of the body, immediately consult your physician without any delay!
- The device contains a Lithium-Ion rechargeable battery that might swell and deform the housing of the device. If you notice any swelling or deformation in the device, stop using it immediately and do not attempt to replace the

swollen battery by yourself. Swollen batteries can be a serious safety hazard, posing risks such as leakage, overheating, or even explosion.

- Do not attempt to remove, replace, charge, puncture or modify the swollen battery. Do not ignore signs of battery damage. Using the device with a damaged battery can affect its performance and permanently damage the device. Avoid handling the device with bare hands if it shows signs of damage or leakage, and do not open the battery compartment if leakage occurred. Avoid direct contact with the leaked electrolyte. If contact occurs, rinse with water and seek medical attention if needed.
- If you have any questions or concerns, seek professional guidance at your hearing care professional before using the device.
- Ship the rechargeable devices according to local regulations as the devices contain lithium-ion batteries. For shipping, they are classified as dangerous goods. If in doubt, ask the courier responsible for shipping how to ship the devices properly.
- Charge the battery only when inserted into the device. Charging the battery otherwise might result in overheating of the battery.
- Dispose of the device or replaced battery at waste disposal sites intended for electronic waste or give your device/battery to your hearing care professional for disposal.
- To supply power to the device, use only IEC 60950-1*, IEC 62368-1 or IEC60601-1 certified charging equipment (e.g. power supply).
- Do not place these rechargeable devices in your checked-in luggage when traveling by plane as they contain lithium-ion batteries. The devices

should be placed in your carry-on/hand luggage.

- Don't charge your device at temperatures outside the operating conditions indicated in the User Guide (see chapter Specifications). Otherwise the charging reliability cannot be guaranteed.
- This device is not for children below 36 months. It contains small parts that can cause choking, if swallowed by children. Keep out of reach of children and mentally challenged people or pets. If swallowed, consult a physician or hospital immediately.
- Keep your device at least 15 cm (6 inches) away from the active implant. If you experience any interference do not use your device and contact the manufacturer of the active implant. Please, note that the interference can also be caused by powerlines, electrostatic discharge, airport metal detectors etc.
- Protect the device connectors, plugs, power supply from dirt and debris. Never immerse them in water.
- Never wash the microphone inputs. Doing so could cause it to lose its special acoustic features.
- Do not place the charging accessories close to an inductive cooking surface. Conductive structures inside the charging accessories could absorb inductive energy resulting in thermal destruction.
- When storing the replaceable battery outside the device, don't expose it to conditions other than indicated in the User Guide (see chapter Information and explanation of symbols) for storage and transportation. This prevents the properties of the battery from deteriorating.

- If the device is not used for a long time, please charge the battery to approximately 50%, remove the battery from the device and keep it at conditions as indicated in the User Guide (see chapter Information and explanation of symbols) for storage and transportation. This mitigates the risk that the battery is deep-discharged by the device. For safety reasons, a deep-discharged battery cannot be re-charged and has to be replaced.
- To replace the battery please follow the instructions provided in the user guide (chapter Care and maintenance).

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* this standard may not be valid anymore in your country, please refer to IEC 62368-1

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029-3384-49/V1.00/2026-02/km © 2026 Sonova AG All rights reserved